



Technical support: support@abbkine.com

Website: https://www.abbkine.com

## Cleaved-Tumstatin (P1426) Polyclonal Antibody

Cat #: ABP50057 Size: 30µl /100µl /200µl

## **Product Information**

	Product Name: Cleaved-Tumstatin (P1426) Polyclonal Antibody		
	Applications: WB, ELISA		Isotype: Rabbit IgG
	Reactivity: Human, Monkey		
REF	Catalog Number: ABP50057	LOT	Lot Number: Refer to product label
	Formulation: Liquid		Concentration: 1 mg/ml
ĵy	Storage: Store at -20°C. Avoid repeated	$\Lambda$	Note: Contain sodium azide.
1	freeze / thaw cycles.	<u>د: ۲</u>	


Background: Type IV collagen, the major structural component of basement membranes, is a multimeric protein composed of 3 alpha subunits. These subunits are encoded by 6 different genes, alpha 1 through alpha 6, each of which can form a triple helix structure with 2 other subunits to form type IV collagen. COL4A3 encodes alpha 3. In the Goodpasture syndrome, autoantibodies bind to the collagen molecules in the basement membranes of alveoli and glomeruli. The epitopes that elicit these autoantibodies are localized largely to the non-collagenous C-terminal domain of the protein. A specific kinase phosphorylates amino acids in this same C-terminal region and the expression of this kinase is upregulated during pathogenesis. COL4A3 is also linked to an autosomal recessive form of Alport syndrome. The mutations contributing to this syndrome are also located within the exons that encode this C-terminal region. Like the other members of the type IV collagen gene family, COL4A3 is organized in a head-to-head conformation with another type IV collagen gene so that each gene pair shares a common promoter.

<u>Application Notes</u>: Optimal working dilutions should be determined experimentally by the investigator. Suggested starting dilutions are as follows: WB (1:500-1:2000), ELISA (1:40000). Not yet tested in other applications.

Storage Buffer: PBS containing 50% Glycerol, 0.5% BSA and 0.02% Sodium Azide.

**Storage Instructions:** Stable for one year at -20°C from date of shipment. For maximum recovery of product, centrifuge the original vial after thawing and prior to removing the cap. Aliquot to avoid repeated freezing and thawing.



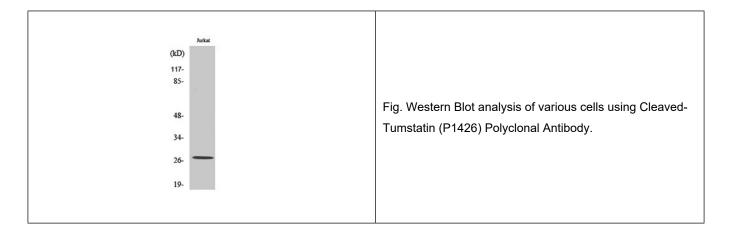

**Note:** The product listed herein is for research use only and is not intended for use in human or clinical diagnosis. Suggested applications of our products are not recommendations to use our products in violation of any patent or as a license. We cannot be responsible for patent infringements or other violations that may occur with the use of this product.

